




Part of the SA3000 Care Call System

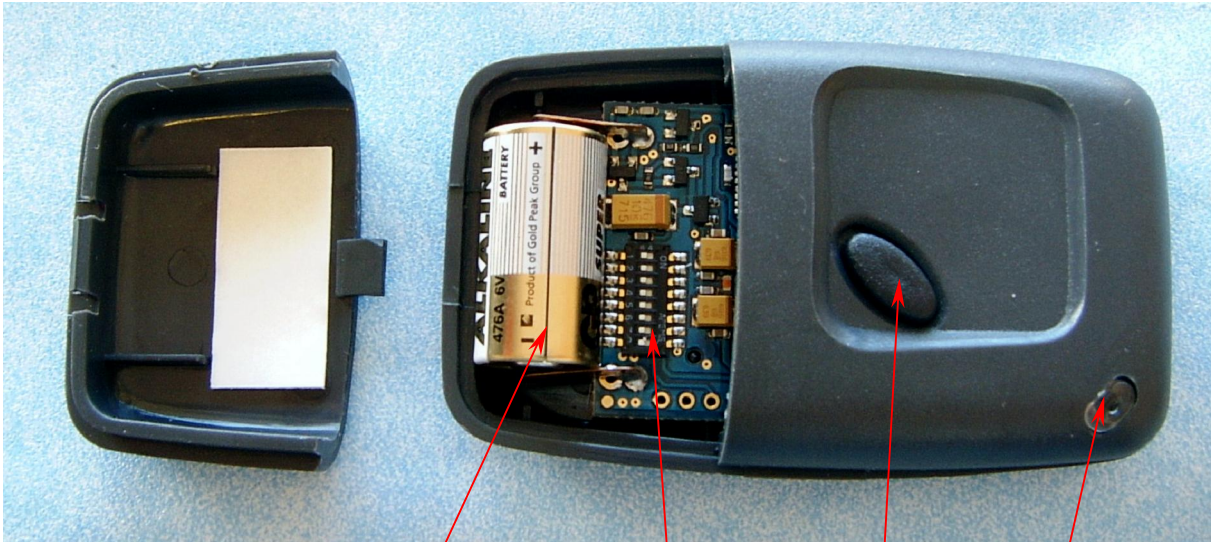

Battery in Battery  
Compartment

Code Switches

Send Button

Light

## Getting Started

1. Ensure the pager unit is fully charged and either alarm clock charger or trickle charger are installed.
2. Remove battery cover using a small screwdriver & install a 6 Volt Silver Oxide (GP476, 4SR44, 4G13, PX28) battery.
3. Press the send button, the green light should shine on the Key Fob and the Pager should vibrate and key 4 will illuminate.
4. This device also incorporates a fall alarm. When the device accelerates rapidly or is shaken roughly it will automatically send a radio signal to the pager.
5. The first 4 of the 8 miniature keys adjacent to the battery are for system coding (On away from the battery). Key 5 will change the zone to 1, key 6 to zone 2 and key 7 to zone 3.
6. A silver Oxide battery should have a large number of activations. The green light should shine for approximately one second. If the battery is low the orange light will shine and a low battery signal will be sent.
7. There is a lockout period of 12 seconds after the send button has been pressed.

**DO NOT IMMERSE THE KEY FOB IN WATER OR ANY OTHER LIQUID. DO NOT SUBJECT THE UNIT TO EXCESSIVE MECHANICAL SHOCK.**

9 100% Recycled Paper

**New 2 Year Warranty**

On All Care Call Products produced after 15th January, 2010

The SA3000 System comprises a Pager Unit, either an Alarm Clock Charger or Trickle Charger, or SignWave and at least one Monitor, Key Fob or Movement Monitor. Clofield Ltd.

Help Line 0800-387-397 Txt & Fax 01246-450789 e-mail [care.info@silent-alert.co.uk](mailto:care.info@silent-alert.co.uk)