


# CCM3A-2216-EU-PIR Care Call Movement Monitor

Part of the SA3000 Care Call System

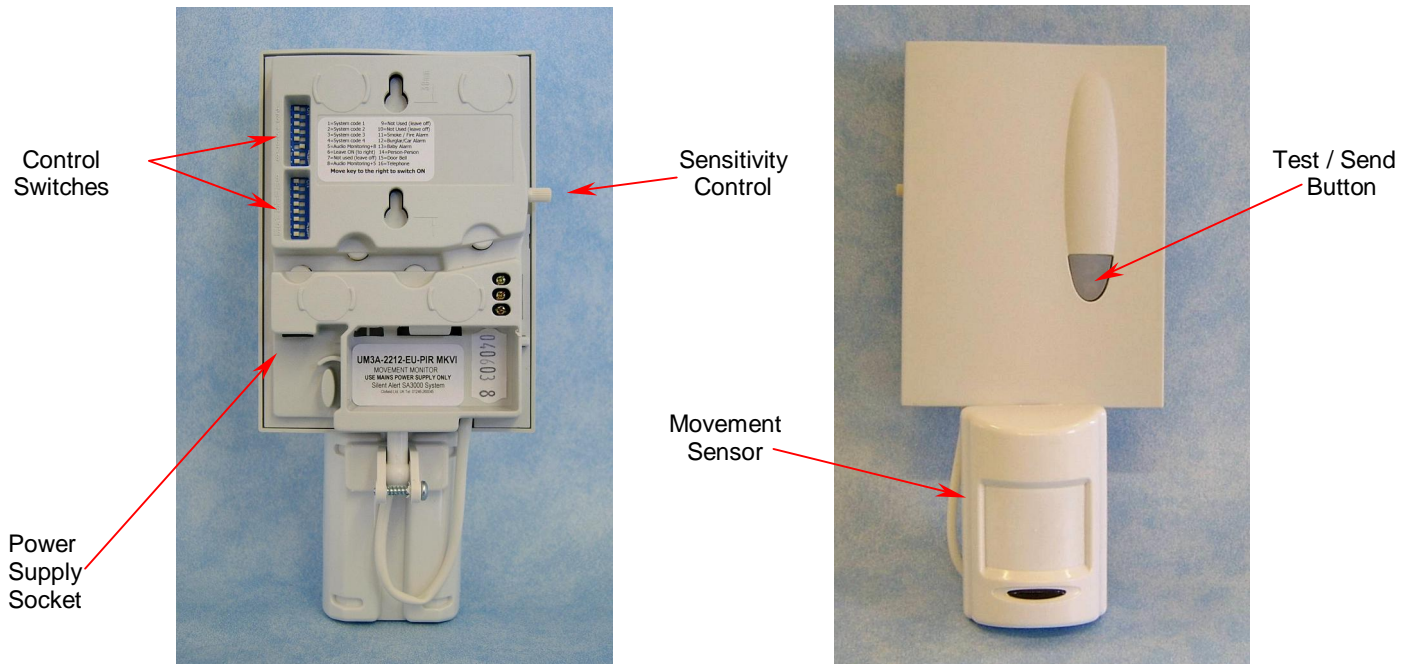

## Getting Started

1. Ensure the Pager unit is fully charged and the Trickle Charger is installed.
2. Plug the Power Supply into a 13amp mains socket.
3. The test button will flash orange until the PIR red light goes out. This is normal.
4. Move one of the control switches (11-16) to the right to determine which key lights on the Pager (See Table).
5. Angle the movement Sensory to the desired position and test the monitor in location before using the fixing screws.
6. Control switches 1-4 are for system coding (On to the right) switch 9 is used to create zones 7 to 12 see table below.
7. To activate the sound monitoring feature move control switch 5 to the On position (to the right) and adjust the sensitivity control to trigger at the desired sound level. It is possible to trigger the sound activation feature by handling the Monitor or whilst adjusting the sensitivity. There is also a 12 second delay after activation before the Monitor can be re triggered. Control switch 7 introduces a short delay before the audio triggers the unit. Switch 8 introduces a medium delay and 7 & 8 together a long delay. An amber light shines when the audio monitor detects sounds. When a sound exceeds a preset delay, the amber light goes out and the green light shines to confirm that the signal has been sent.

Key on Pager	Control Switch
1	16
2	15
3	14
4	13
5	12
6	11

Key on Pager	Control Switch
7 (Battery+1)	9+16
8 (Battery+2)	9+15
9 (Battery+3)	9+14
10 (Battery+4)	9+13
11 (Battery+5)	9+12
12 (Battery+6)	9+11

**New 2 Year Warranty**  
On All Care Call Products produced after 15th January, 2010

9 100% Recycled Paper

**DO NOT IMMERSER THE UNIVERSAL MONITOR IN WATER OR ANY OTHER LIQUID. DO NOT SUBJECT THE UNIT TO EXCESSIVE MECHANICAL SHOCK.**