

# Improved Silent Alert Care Call System with Pager and SignWave receiver options



Wireless Alerting System for Carers



# The Silent Alert Care Call System is a radio paging system used by carers at home or in a care facility



Ideal for the care of people with

- Alzheimer's
- Mobility Problems
- A tendency to fall
- A tendency to wander or sleepwalk

# Vibrating Pager Option



A small Pager unit worn by the Carer. When the Pager receives a signal it illuminates one of the numbered keys and gently vibrates to discreetly alert the carer of an event

The Pager can respond to a maximum of 12 transmitters allowing one Carer to monitor multiple patients

If a signal is received at night, whilst the pager is charging, an audio alarm or vibrating pad will wake the Carer



# SignWave™ Portable Sounder Flasher option



6

5

4

3

2

1



Incorporates a miniature high intensity strobe & ultra loud sounder with adjustable volume

Standard open air operating range of 100m, with a long range setting of 500m at the flick of a switch



6

5

4

3

2

1

The Person to Person Key Fob can be used to page a Carer who is in another part of the home or care facility

Ideal when caring for someone who has mobility problems or is prone to falling



Now with shock detection option

# Mini Magnetic Door Monitor



6

5

4

3

2

1



The Mini Magnetic Door Monitor transmits when a door is opened

Ideal for alerting the carer when and where a patient moves around the house



# Chair, Bed and Carpet Pressure Pads



6

5

4

3

2

1



Ideal for use when Carers discreetly wish to monitor a patient leaving their bed or chair



# Movement Monitor for someone who sleepwalks or has Alzheimer's



6

5

4

3

2

1



Locate on the landing  
or at a bedside

Can be triggered by  
either movement,  
sound or both

New optional delay prevents  
bathroom visits triggering the unit

If the patient does not return  
within a pre-set time the unit will  
alert the carer



# Audio Monitor



6

5

4

3

2

1



The Universal Monitor allows the user to monitor a room or area by responding to predetermined levels of sound

When a sound of the correct volume and duration is picked up by the monitor, it will send a signal to the Pager or Signwave



# Special Monitors




Our Special monitors are adapted to interface with popular third party products to enable them to be monitored with our Pager or SignWave



Can be linked to:

- EMFIT Epilepsy Monitor
- Bed occupancy monitor
- Neonatal Respiration Monitor

Dual range facility of 100m, adequate for most homes and gardens and 500m for use within a care facility or when caring for a neighbour

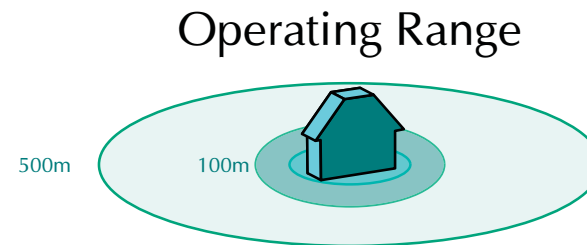

Up to 12 zone options using a multiplier key