

# Introducing the Silent Alert Care Call System



Wireless alerting system for Carers.  
Ideal for stand alone telecare and personal  
alarm use.



# Care Call is a radio paging system used by carers at home or in a care facility



With DIY installation and setup the Care Call system fills the gap between entry level border alarms and expensive nurse call station systems. The ideal solution for caring for 1-12 people per receiver.

Ideal for the care of people with

- Alzheimer's
- Mobility problems
- A tendency to fall
- A tendency to wander or sleepwalk



# Care Call Pager




When the Pager receives a signal it illuminates one of the numbered keys and beeps or discreetly vibrates to alert the Carer of an event. The units internal battery is recharged each night in the Care Call charger.

The Pager can respond to 12 transmitters allowing one Carer to monitor multiple areas or patients.

If a signal is received at night, whilst the pager is charging, an audio alarm or vibrating pad will wake the Carer.



Operating Range

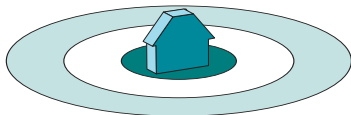

The Care Call System has a range of up to 1000m in open air



Made in the UK

# SignWave™ portable sound & flash receiver



- Easy to use cancel button
- LED Strobe light
- Volume Control
- Back lit display
- Up to 1000M operating range
- Simple 'drop in' charging base

The SignWave incorporates a high intensity LED strobe & loud (100dB) sounder with adjustable volume control. The unit can be placed anywhere in the home to alert Carers to an event.

Ergonomically designed, the unit is easy to pick up and operate with its large cancel button on the top of the unit.

Discreet overnight cover can be achieved using a vibrating pillow pad attached to the charging base.



Operating Range



The Care Call System has a range of up to 1000m in open air



Made in the UK

# Care Call

## Mini Magnetic Door Monitor



The Mini Magnetic Door Monitor transmits when a door is opened. Ideal for alerting the Carer when and where a patient moves around the house.

Up to 12 doors can be monitored by using additional monitors.

The door monitor can be configured with a trigger delay of 10, 20 or 30 seconds.



# Chair, bed leaving and under carpet pressure pads

- 1
- 2
- 3
- 4
- 5
- 6
- +6

Ideal for use when Carers wish to discreetly monitor a patient leaving their bed or chair.



# The Person to Person key fob can be used to page a Carer who is in another part of the home or care facility



Ideal when caring for someone who has mobility problems or is prone to falling.

The keyfob can be worn around the neck. A simple press of the button will alert the carer.

The G-FOB has an internal sensor that can help detect falls. When a shock is detected the keyfob will automatically transmit to alert the Carer.



# Movement Monitor for someone who is prone to wandering



Simply locate the monitor on the landing or at a bedside.

The movement monitor will be triggered by either movement, sound or both.

New optional return delays prevent bathroom visits from triggering the unit. If the patient does not return within a pre-set time the unit will alert the Carer.





# Care Call Universal Monitor with audio detection



The Universal Monitor allows the user to monitor a room or area by responding to predetermined levels of sound.

When a sound of the correct volume and duration is picked up by the monitor, it will send a signal to the Pager or Signwave.



# Care Call smoke and carbon monoxide alarms

- 1
- 2
- 3
- 4
- 5
- 6
- +6



Smoke and CO alarms can be added to any system to alert the Carer to an emergency. Ionisation and optical smoke alarms are available.



# Care Call special monitors




Our Special monitors are adapted to interface with popular third party products to enable them to be monitored with our Pager or SignWave  
Can be linked to:


- EMFIT Epilepsy Monitor
- Bed occupancy monitor
- Neonatal Respiration monitor



Dual range facility of 100m, adequate for most homes and gardens and up to 1000m for use within a care facility or when caring for a neighbour.

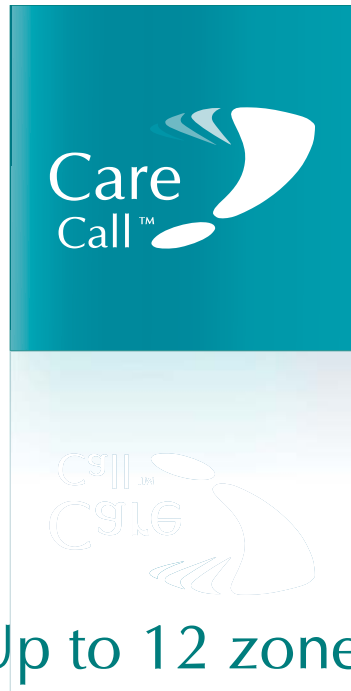

Operating Range

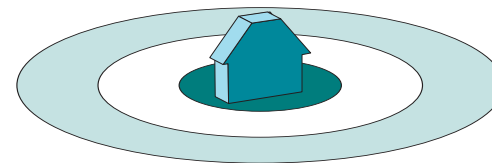

The Care Call System has a range of up to 1000m in open air

Up to 12 zone options using a multiplier key

# Systems priced from just £175.00

•  
Free telephone support

•  
1 Year Warranty

